


Proscale Orion 150 Pocket Scale

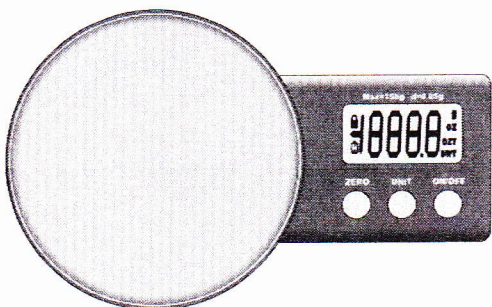

CARE AND MAINTENANCE

1. Do not dismantle the scale.
2. The scale should not be used while vibrating or in direct sunlight.
3. Rapid change of temperature will affect the accuracy of the scale. However, the scale will compensate about two hours after the change of temperature.
4. Do not spill liquid on the scale. You may wipe the scale down periodically with a damp cloth.
5. Keep the scale away from heat source, ovens, stoves, etc.
6. Keep the scale in a cool, dry, dustless environment and remove batteries in case of extended storage.
7. Replace batteries in scale before returning your scale for service. If you find service is still required, remove the batteries and repack the scale in its original carton and packing materials for transportation.

This scale is a precision instrument and specially adjusted and manufactured for correct measuring. You are requested to use the scale as per the following instructions.

USE

1. Place the scale on a hard, level surface.
2. Push ON/OFF button. The scale will display **-18888**, for about three seconds then it will display **0.00**. The scale is now ready to weigh.
3. Place the object(s) to be weighed on the scale. The display will show the weight. Should the item(s) being weighed surpass the scale capacity of 150g, **0-Ld** will appear on the display. This means over load error.
4. Zero system : Activate the scale by pressing the ON/OFF button, press the ON/OFF button again and the display will read **0.00**. Place the first item to be weighed individually by repressing the ON/OFF button to **0.00** each time an item is added until the 150g capacity is reached.
5. UNIT system: Press the UNIT button will change the unit of the scale.
6. Power off: After load is stabilized and measuring is completed, scale shuts itself off in thirty seconds. (prolonging battery life) Press the ON/OFF buttons will shut it off too.
7. Battery replacement: When batteries need to be replaced (**LO**), will appear on the display, Press the off button and replace the batteries immediately.

CAUTION

The scale may not display 0.00 if a weight of more than 150 grams is placed on the tray or if the scale receives a hard shock / drop. If this happens please try pressing the on/off key to reset the scale to 0.00. You may also try recalibration. If the scale is not going to be used for an extended period of time please remove the batteries to avoid possible leakage.

Calibration instructions

1. Place the scale on a flat, very stable surface and power it OFF.
2. Press "ON/OFF" key, turning on the scale. The LCD will display "0.00".
3. Press and hold the "UNIT" key, the display will show "CAL". Release the key.
4. Press "UNIT" key again, the LCD will display "100.00", gently put a 100g weight on the platform. The display will show "PASS".
5. Calibration is complete. Remove the weight and turn the scale off.

www.proscale.net